



A tissue processing **revolution**


Revos Tissue Processor from Epredia

Syntri Safeguard Consumables from Epredia



With the Revos Rotational Tissue Processor,

your workflow can experience
multiple improvements and provide
your laboratory confidence with high
quality tissue.

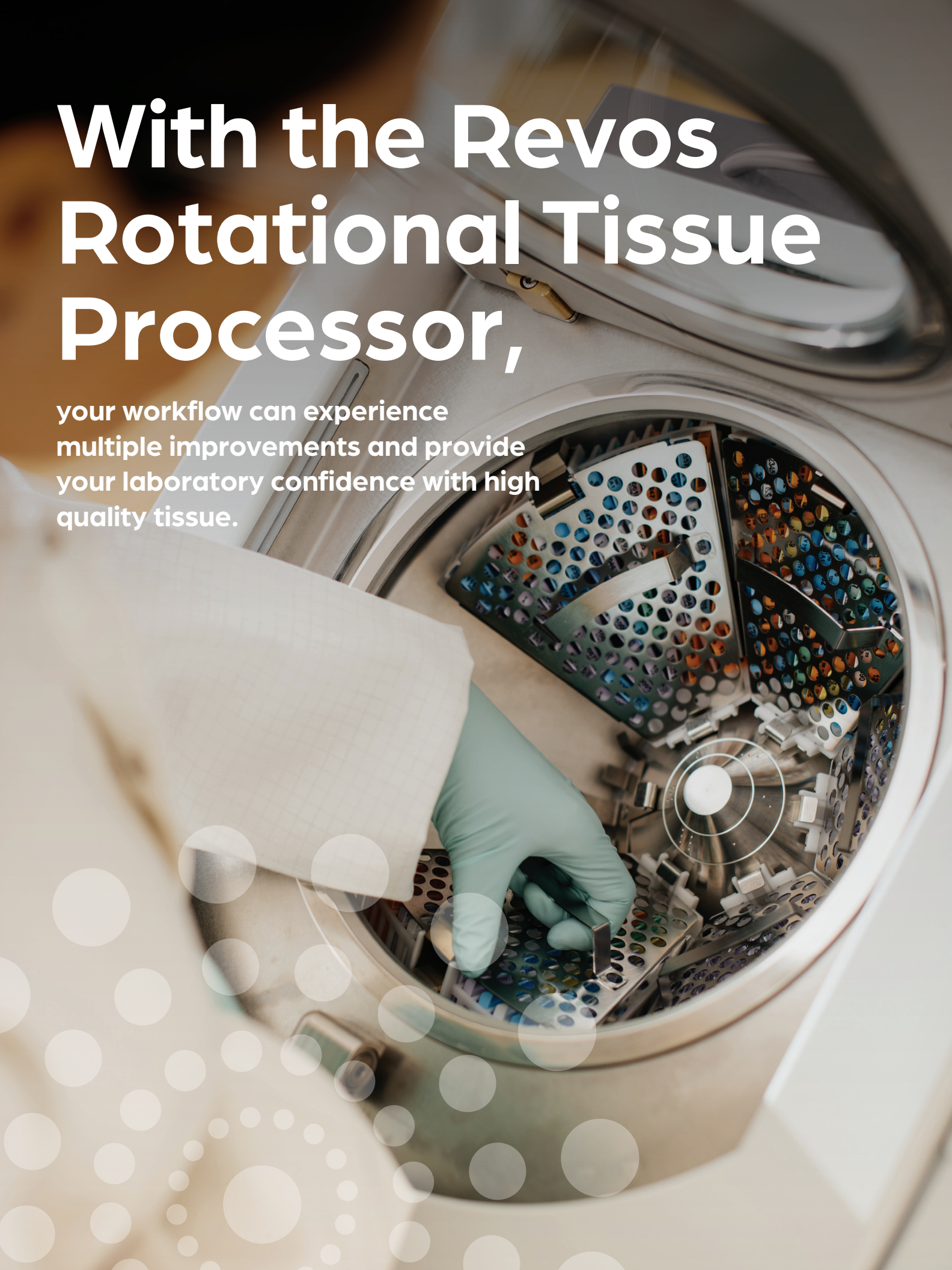

Revos tissue processor

Bottlenecks in laboratory workflows reduce throughput, lengthen the time from tissue sampling to diagnoses and increase costs.

The Revos™ workflow-enhancing tissue processor, paired with Syntri Safeguard™ consumables, overcomes these challenges. It offers reduced tissue processing time and less risk of tissue damage – all with no reduction in processing quality.



Rapid processing

Our unique, canted chamber enhances the distribution of reagent, reduces tissue processing time and allows for high-quality processing results.



Intelligent monitoring



The instrument's intelligent reagent management system tracks and logs reagent usage. It remembers every bottle loaded, protocol executed and service performed – allowing laboratories to use reagents more efficiently.

Safe environment

The Revos instrument's dual-filtration system and down-draft ventilation helps protect laboratory staff by controlling both solvent and formalin fumes.

Protecting tissue

Be ready for unexpected power outages – the onboard battery backup system keeps your samples safe and protocols active for up to three hours.

Syntri enabled

Syntri Safeguard technology enables the instrument to sense where the 100% processing alcohol is loaded in the system and prevents incorrect replenishment

VNC Remote

Revos Tissue Processor allows the laboratories' user interface to be viewed and interacted with remotely via a virtual connection. This remote access to the instrument allows our certified, trained field service technicians to quickly identify the problem and promptly get you back to uptime with VNC (Virtual Network Computing) Remote.

Laboratories are under increasing pressure to speed up processes, reduce costs and eliminate errors. The Revos tissue processor and Syntri Safeguard consumables are expressly designed to work together and help overcome these challenges.





- 1 Color touch screen monitor
- 2 300-cassette-capacity chamber
- 3 Reagent storage cabinet
- 4 Downdraft ventilation
- 5 Carbon and potassium permanganate filtration
- 6 USB port for easy downloading of data

Revos workflow-enhancing tissue processor

Specifications

Dimensions (D x W x H)	580 mm x 850 mm x 1370 mm 23 in x 33.5 in x 54 in
Weight:	• Instrument = 165 kg / 364 lbs • Instrument + typical reagents = 240 kg / 530 lbs
Power requirements	• 100-240 V AC • 1300 W / 110 V AC = 11.8 A • 1300 W / 230 V AC = 5.6 A
Backup battery	Up to 3 hours
Certifications	• IEC61010 3rd Ed. • UL • CE
Number of retorts	1
Cassette capacity	• Standard 300 • Super Mega 42
Alcohol configuration	6 bottles
Xylene/Xylene Substitute configuration	3 bottles
Fixation configuration	2 bottles
Number of protocols	10
Data connection	LIMS connectivity via RS232 port

Ordering Information

REF Number

Instrument

Revos Tissue Processor

A84100001ATC
A84100001

Filters

Formaldehyde (potassium permanganate) filter; 6 per pack

B9990612CS

Charcoal filters; 6 per pack

7411258

Other accessories

Downdraft extraction duct adapter kit

A84110045

Main air system duct adapter kit

A78410024

Storage bottles (5L); 6 per pack

A78410026

Waste wax drawer; 5 per pack

8300



A smarter solution

The accuracy of patient tests often hinges upon proper preparation and consistency in processes throughout the laboratory.

When paired with Revos™ reagents featuring Syntri Safeguard™ RFID technology, the Revos instrument tracks and logs reagent usage. It remembers every bottle loaded, protocol executed and service performed – helping to reduce downtime and facilitating efficient use of reagents.

Syntri Safeguard Consumables

Description	Cat. Number
Syntri Safeguard™ 10% Neutral Buffered Formalin, 4 x 5 L	5735SS
Syntri Safeguard™ Pen-Fix™, 4 x 5 L	6115SS
Syntri Safeguard Dehydrant 100%, 4 x 5 L	6215SS
Syntri Safeguard Xylene, 4 x 5 L	6615SS
Syntri Safeguard™ Clear-Rite™ 3, 4 x 5 L	6915SS
Syntri Safeguard™ Flex™ 100, 4 x 5 L	8115SS
Syntri Safeguard Isopropyl Alcohol, 4 x 5 L	9517SS
Syntri Safeguard™ Formal-Fixx™, 4 x 5 L	6764215SS
Shandon™ Syntri Safeguard™ Xylene Substitute, 4 x 5 L	6764515SS

Syntri Safeguard Consumables

4481 Campus Drive
Kalamazoo, MI 49008
UNITED STATES
+1 (800) 522-7270

Revos Tissue Processor and Equipment

 Tudor Road, Manor Park
Runcorn, WA7 1TA
UNITED KINGDOM
+44 (0) 800 018 9396
+44 (0) 1928 534 000

Find out more at [epredia.com](https://www.epredia.com)

 Enhancing precision
cancer diagnostics